




DEBUNK Plus+ Fluorescence Immunoassay Analyzer AFR-700

- High Accuracy (Stringent Quality Control Monitoring & System Check Function)
- User-Friendly (Simple Plug & Play Operation, 10.1" LCD Touch Screen)
- Easy Operation (Less Training Required)
- High Efficiency (Batch Testing, 12 Tests in One Go)
- Time-saving (Reading Time less than 20 Seconds)
- Get Printing Test Reports Instantly (Built-in Thermal Printer)
- LIS/HIS Compatibility Available
- Strong Stability
- Extensive Test Menu
- Clear Status Indicators (5 different color indicators for easily monitoring of the testing process)
- Effortless Standard Curve Import (Uses RFID cards for the automatic induction)
- Automated Recognition (Automatically identifies QR codes on cassettes)
- Incubation Temperature Control (Ensuring precise incubation at 37°C)

Product Features:

Powered by Europium-based Technology for High Quality
Fast Result (3-15 minutes)
Simple Operation (Less Training Required)
Objective (Results Read by Analyzer)
Stringent Quality Control Insure High Accuracy
Easy Storage Condition (4-30°C) with 2 Years Shelf Life

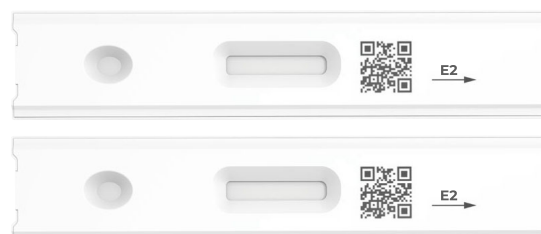


Test Menu

Cat. No.	Product	Specimen	Cut-off	Measuring Range	Pack	CE status
Inflammation Marker						
FI-CRP-402	CRP	Fingerstick WB, WB/S/P	hsCRP: < 1 mg/L, CRP: < 10 mg/L	0.5-200 mg/L	10T/25T	CE
FI-PCT-402	PCT	Fingerstick WB, WB/S/P	< 0.5 ng/mL	0.1-50 ng/mL	10T/25T	CE
FI-IL6-402	IL-6	WB/S/P	< 10 pg/mL	3.0-5000 pg/mL	10T/25T	CE
FI-SAA-402	SAA	Fingerstick WB, WB/S/P	< 10 mg/L	1-300 mg/L	10T/25T	CE
Tumor Marker						
FI-PSA-302	PSA	S/P	< 4 ng/mL	2-100 ng/mL	10T/25T	Non-CE
FI-PSA-402	PSA	Fingerstick WB, WB/S/P	< 4 ng/mL	0.1-100 ng/mL	10T/25T	Non-CE
FI-FPSA-302	Free PSA	S/P	/	0.1-30 ng/mL	10T/25T	Non-CE
FI-AFP-302	AFP	S/P	< 20 ng/mL	5-400 ng/mL	10T/25T	CE
FI-CEA-302	CEA	S/P	< 4 ng/mL	1-500 ng/mL	10T/25T	CE
FI-FOB-602	FOB	Feces	50 ng/mL	10-1000 ng/mL	10T/25T	CE
FI-CA211-402	CYFRA21-1	Fingerstick WB, WB/S/P	<3.3 ng/mL	1-300ng/mL	10T/25T	CE
Cardiac Marker						
FI-CTI-402	cTnI	Fingerstick WB, WB/S/P	< 0.5 ng/mL	0.1-40 ng/mL	10T/25T	CE
FI-CTNT-402	cTnT	Fingerstick WB, WB/S/P	< 0.5 ng/mL	0.2-40 ng/mL	10T/25T	CE
FI-NBNP-402	NT-proBNP	Fingerstick WB, WB/S/P	< 0.45 ng/mL	0.3-22 ng/mL	10T/25T	CE
FI-CKMB-402	CK-MB	Fingerstick WB, WB/S/P	< 5 ng/mL	0.2-75 ng/mL	10T/25T	CE
FI-HFA-402	H-FABP	Fingerstick WB, WB/S/P	< 8 ng/mL	1-120 ng/mL	10T/25T	CE
FI-MYO-402	Myoglobin	Fingerstick WB, WB/S/P	< 90 ng/mL	5-200 ng/mL	10T/25T	CE
FI-AMI-435	Troponin I/Myoglobin /CK-MB (3 in 1)	Fingerstick WB, WB/S/P	cTnI: < 0.5 ng/mL MYO: < 90 ng/mL CK-MB: < 5 ng/mL	0.1-40 ng/mL 5-200 ng/mL 0.2-75 ng/mL	10T/25T	CE
FI-AMIT-435	Troponin T / Myoglobin / CK-MB Combo	Fingerstick WB, WB/S/P	cTnT: < 0.2 ng/mL MYO: < 90 ng/mL CK-MB: < 5 ng/mL	0.1-40 ng/mL 5-200 ng/mL 0.2-75 ng/mL	10T/25T	CE
FI-CPC-435	Troponin I/D-Dimer/ NT-proBNP (3 in 1)	Fingerstick WB, WB/S/P	cTnI: < 0.5 ng/mL D-Dimer: < 0.5mg/L NT-proBNP: < 0.45 ng/mL	0.1-40 ng/mL 0.1-10mg/L 0.3-22 ng/mL	10T/25T	CE
Coagulation						
FI-DDM-402	D-Dimer	Fingerstick WB, WB/P	0.5 mg/L	0.1-10 mg/L	10T/25T	CE
Endocrinology						
FI-BHCG-402	β-hCG	Fingerstick WB, WB/S/P	< 5 mIU/mL	2-200,000 mIU/mL	10T/25T	CE
FI-BHCG-202	β-hCG	Fingerstick WB, WB/S/P/Urine	< 5 mIU/mL	5-100,000 mIU/mL	10T/25T	CE
FI-FSH-402	FSH	Fingerstick WB, WB/S/P	/	1-100 mIU/mL	10T/25T	CE
FI-LH-402	LH	Fingerstick WB, WB/S/P	/	1-300 mIU/mL	10T/25T	CE
FI-AMH-402	AMH	WB/S/P	/	0.05-25 ng/mL	10T/25T	CE
FI-TES-302	Testosterone	S/P	Female:0-2 ng/mL, Male:2.6-12 ng/mL	1.0-20 ng/mL	10T/25T	CE
FI-P4-302	Progesterone (P4)	S/P	/	1.4-60 ng/mL	10T/25T	CE
FI-E2-402	Estradiol (E2)	Fingerstick WB, WB/S/P	Male:11-44 Female(Follicular period):31~91 Female(Ovulation):60-533 Female(Luteal period):60-233 Female(Pregnant):>153 Female(Menopause):<140	10-3000 pg/mL	10T/25T	CE
FI-PRL-402	Prolactin (PRL)	Fingerstick WB, WB/S/P	Female (Non-pregnant): 2- 25 ng/mL Female (Pregnant): 10 - 209 ng/mL Male:2 - 18 ng/mL	1.0 - 200 ng/mL	10T/25T	CE
FI-T3-402	T3	Fingerstick WB, WB/S/P	1.23-3.08nmol/L	0.62-9.24 nmol/L	10T/25T	CE
FI-FT3-402	FT3	Fingerstick WB, WB/S/P	2.6-5.7 pmol/L	1.5-46 pmol/L	10T/25T	CE
FI-T4-402	T4	Fingerstick WB, WB/S/P	65-155nmol/L	12.87-300 nmol/L	10T/25T	CE
FI-FT4-402	FT4	Fingerstick WB, WB/S/P	9-19.5pmol/L	1-100 pmo/L.	10T/25T	CE
FI-TSH-402	TSH	Fingerstick WB, WB/S/P	0.27-4.2 µIU/mL	0.1-200 µIU/mL	10T/25T	CE
FI-FER-402	Ferritin	Fingerstick WB, WB/S/P	Male:30-400 ng/mL, Female:13-150 ng/mL	1-2000 ng/mL	10T/25T	CE
FI-COR-302	Cortisol	S/P	140 - 690 nmol/L	50 nmol/L - 1000 nmol/L	10T/25T	CE
Nephrology						
FI-ALB-102	mAlb	Urine	< 20 mg/L	5-300 mg/L	10T/25T	CE
FI-CYSC-402	CysC	Fingerstick WB, WB/S/P	< 1.1 mg/L	0.1-9.0 mg/L	10T/25T	CE
FI-BMG-402	β2MG	Fingerstick WB, WB/S/P	< 2.2 mg/L	0.2-20 mg/L	10T/25T	CE
FI-NGAL-102	N-GAL	Urine	< 10 ng/mL	10-1500 ng/mL	10T/25T	CE
Diabetes						
FI-HBA-402	HbA1c		Fingerstick WB, WB/S/P	4-6% 4-14.5%	10T/25T	CE

Cat. No.	Product	Specimen	Cut-off	Measuring Range	Pack	CE status
Gastrointestinal Diseases						
FI-HP-602	<i>H. pylori</i> Antigen	Feces	< 2.5 ng/mL	2-1500 ng/mL	10T/25T	CE
FI-CAL-602	Calprotectin	Feces	20 µg/g	1-1000 µg/g	10T/25T	CE
FI-AD-602	Adenovirus	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-ROT-602	Rotavirus	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-RAC-602	Rota/Adeno	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-NOV-602	Norovirus GI/GII	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-CDG-602	<i>Clostridium difficile</i> GDH	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-CDT-625	<i>C. difficile</i> Toxin A+Toxin B	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-CD-602	<i>C. difficile</i> GDH/Toxin A/Toxin B	Feces	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
Respiratory Diseases						
FI-STA-502	<i>Strep A</i>	Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-MP-502	<i>Mycoplasma pneumoniae</i> Antigen	Throat Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-SP-102	<i>Streptococcus pneumoniae</i>	Urine	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-LEG-102	<i>Legionella pneumophila</i> Antigen	Urine	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-FLU-502	FLU A+B	Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-RSV-502	RSV	Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-HMPV-502	Human Metapneumovirus (HMPV)	Nasopharyngeal Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-HPIV-502	Human Parainfluenza Virus (HPIV)	Nasopharyngeal Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-NCP-502	COVID-19 Antigen	Swab	< 1 COI	Qualitative (1-1000 COI)	20T	CE
FI-NCP-402	COVID-19 IgG/IgM	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-20 COI)	10T/25T	CE
FI-CSG-402	COVID-19 S-RBD IgG Antibody	Fingerstick WB, WB/S/P	< 1 COI	Semi-quantitative (0.5-100 COI)	10T/25T	CE
FI-CSG-402B	COVID-19 S-RBD IgG Quantitative	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0.5-100 COI)	10T/25T	CE
FI-AD-502	<i>Adenovirus pneumoniae</i>	Nasopharyngeal Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
FI-AR-502	Adeno/RSV	Nasopharyngeal Swab	< 1 COI	Qualitative (0-100 COI)	10T/20T	CE
Vector Borne Diseases						
FI-DEN-402	Dengue IgG/IgM	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-DES-402	Dengue NS1	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
FI-MAL-435	Malaria P.f./P.v./Pan Antigen	Fingerstick WB, WB	< 1 COI	Qualitative (0-20 COI)	10T/25T	CE
FI-TYP-402	Typhoid Antibody	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (1-20 COI)	10T/25T	CE
Hepatitis						
FI-HBSB-402	HBsAb	WB/S/P	< 10mIU/mL	5-1000 mIU/mL	10T/25T	Non-CE
FI-HBSG-402	HBsAg	WB/S/P	<0.5 IU/mL	0.5-200 IU/mL	10T/25T	Non-CE
FI-HCV-402	HCV Antibody	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-100 COI)	10T/25T	Non-CE
FI-HIV-402	HIV 1/2 Antibody	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (1-20 COI)	10T/25T	Non-CE
FI-HIG-402	HIV p24 Antigen	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-100 COI)	10T/25T	Non-CE
FI-SYP-402	Syphilis Antibody	Fingerstick WB, WB/S/P	< 1 COI	Qualitative (0-100 COI)	10T/25T	CE
Other Infectious Diseases						
FI-TET-402	Tetanus Antibody	Fingerstick WB, WB/S/P	< 0.1 IU/mL	0.05-10 IU/mL	10T/25T	CE
FI-ASO-402	Anti-streptolysin O (ASO)	Fingerstick WB, WB/S/P	Adults: < 200 IU/mL Children: <250 IU/mL	25-800 IU/mL	10T/25T	CE
FI-MXG-402	Monkeypox Virus Antigen	Fingerstick WB, WB/S/P, Rash Exudate Swab, Throat Swab	< 1 COI	Qualitative (0-100 COI)	20T	CE
Autoimmunity						
FI-GE-402	IgE	Fingerstick WB, WB/S/P	< 100 IU/mL	20-1600 IU/mL	10T/25T	CE
FI-GA-402	IgA	Fingerstick WB, WB/S/P	0.7-4.3 g/L	0.1-7.0g/L	10T/25T	CE
FI-RFM-402	Rheumatoid Factor(RF) IgM	Fingerstick WB, WB/S/P	< 20IU/mL	8-200 IU/mL	10T/25T	CE
Others						
FI-VD-402	Vitamin D	Fingerstick WB, WB/S/P	>30 ng/mL	5-100 ng/mL	10T/25T	CE
Drugs of Abuse						
FI-AMP-102	AMP	Urine	< 1,000/500 ng/mL	300-3000 ng/mL	10T/25T	CE
FI-BAR-102	BAR	Urine	< 300 ng/mL	50-1200 ng/mL	10T/25T	CE
FI-BZO-102	BZO	Urine	< 300 ng/mL	150-1500 ng/mL	10T/25T	CE
FI-COC-102	COC	Urine	< 300 ng/mL	50-3000 ng/mL	10T/25T	CE
FI-MET-102	MET	Urine	< 1000 ng/mL	100-5000 ng/mL	10T/25T	CE
FI-MTD-102	MTD	Urine	< 300 ng/mL	150-1500 ng/mL	10T/25T	CE
FI-OP1-102	OPI	Urine	< 2000 ng/mL	300-6000 ng/mL	10T/25T	CE
FI-THC-102	THC	Urine	< 50 ng/mL	25-150 ng/mL	10T/25T	CE
FI-APVP-102	α-PVP	Urine	< 300 ng/mL	150-3000 ng/mL	10T/25T	CE
FI-MDMA-102	MDMA	Urine	< 500 ng/mL	250-1500 ng/mL	10T/25T	CE
FI-MOP-102	MOP	Urine	< 300 ng/mL	150-6000ng/mL	10T/25T	CE
FI-PCP-102	PCP	Urine	< 25 ng/mL	10-150ng/mL	10T/25T	CE

Products marked in red : Fingerstick Whole Blood Specimen Available

Fluorescence Immunoassay Analyzer AFR-700

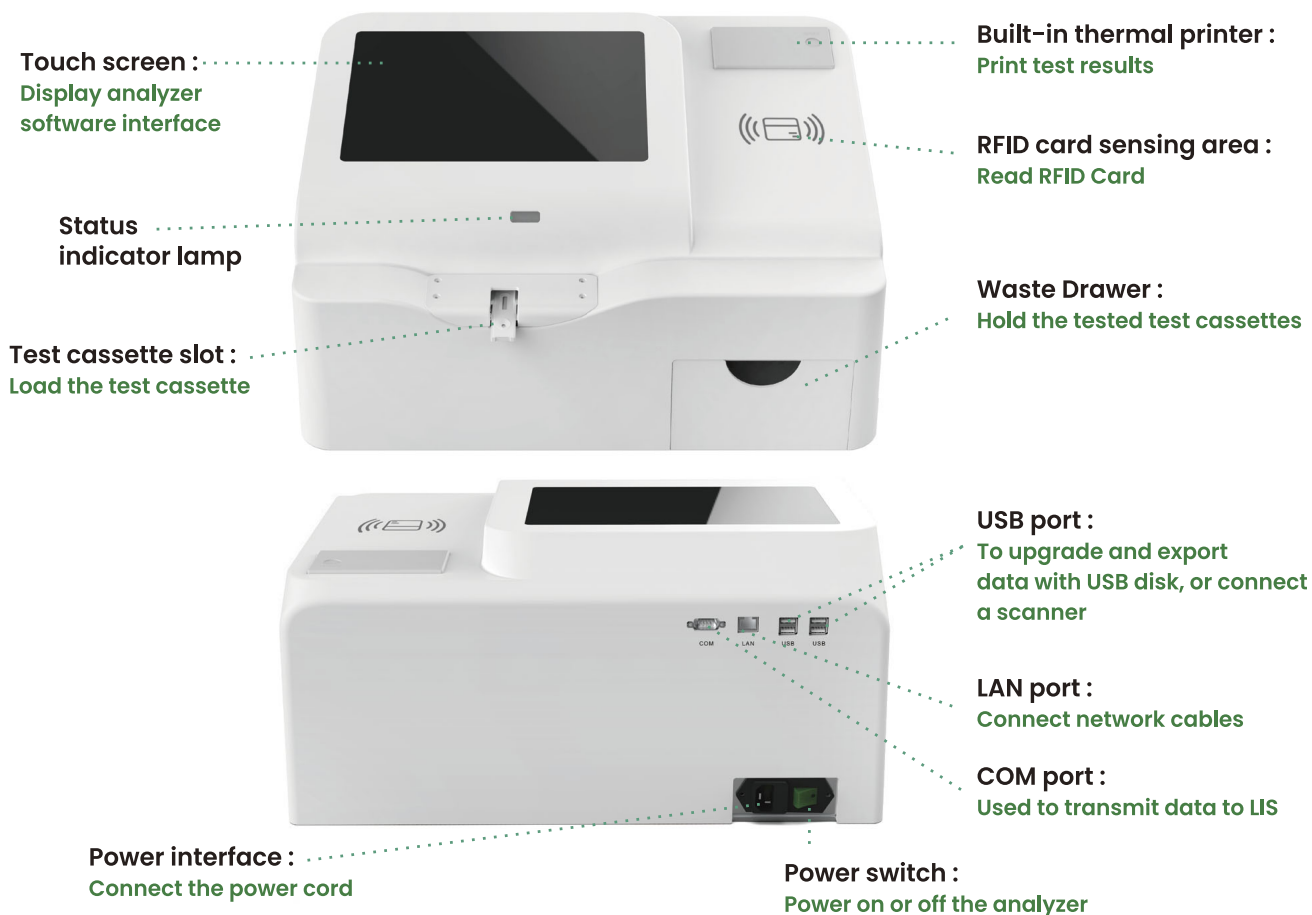

Technical Specifications

Principle	Fluorescence Immunoassay	Working Voltage	AC 100-240V 50/60 Hz
Test Formats	Cassette	Connecting Ports	COM*1 USB *4 LAN *1
Measurement	Quantitative, Qualitative	Memory	10000 records
Read Time	< 20 seconds	Operating Conditions	10-30 °C
Dimensions (L*W*H)	408 mm * 502 mm * 260 mm	Relative Humidity	RH ≤ 85%
Display Dimensions	10.1 inch color LCD	Excitation / Emission band	365/610 (nm)
Weight	13.7 kg	Printer	Built-in Thermal printer
Maximum Input Power	200VA	test channels	12

Distribuidor exclusivo en España

MONLAB

Monlab S.L.

Carrer Cobalto 74, 08940 Cornellà del Llobregat
Barcelona - Spain

Tel. +34 93 433 58 60 - Fax +34 93 436 38 94
mn.sales@monlab.es- www.monlab.com

Hangzhou AllTest Biotech Co.,Ltd.

+86-571-56267891

www.alltests.com.cn

info@alltests.com.cn

No. 550, Yin Hai Street, Hangzhou Economic and Technological Development Area, Hangzhou P.R. China